


Choice of 6 finishes in 8 profiles




Bleached Oak also available for Standard profiles

Standard Melamine Profiles

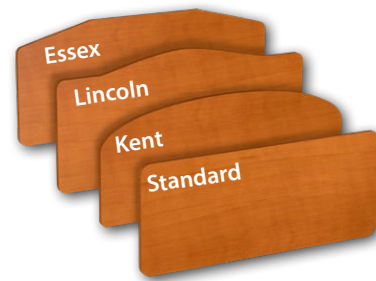

Elite HPL Profiles

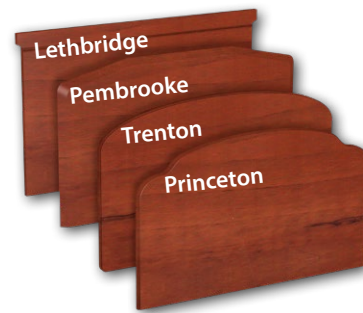

Standard features include:

- Durable, anti-microbial polyester powder-coating of metal components
- Smooth deck design for easier cleaning
- Mobility at any height
- Foot end controller with lockout
- Vascular foot rest with easy ratchet operation
- Beds move vertically with virtually no horizontal motion
- Pendant bed controller

Optional features:

- Advanced mobility as well as floor level locking
- Molded half-head assist rail
- Bed carrier for streamlined bed logistics and maintenance
- Stylish head and foot board options
- Bed extenders



ADVANTAGE
READYWIDE SERIES

Height Floor to Mattress Support	7-7/8" - 31" / 200 mm - 787mm
Deck Angles - Back Elevation	0 - 70°
Deck Angles - Knee Elevation	0 - 30°
Bed Length - Overall Length	82", 86", 90" / 2083 mm, 2184 mm, 2286 mm
Bed Length - Sleep Surface	76", 80", 84" / 1930 mm, 2032 mm, 2134 mm
Bed Width - Sleep Surface	35", 39", 42" / 889 mm, 991 mm, 1067 mm
Weight of Bed (without accessories)	259 lbs / 118 kg (Standard model) 270 lbs / 123 kg (Plus model)
Safe Working Load	500 lbs / 227 kg (Standard model) 650 lbs / 295 kg (Plus model)
Warranty	Mechanical: 5 years Frame: 10 years Welds: 15 years

Item #	Description
MC7W	Advantage ReadyWide with foot end controls
MC7WM	Advantage ReadyWide with mobility & foot end controls
MC7WL	Advantage ReadyWide with floor lock & foot end controls
MC7WML	Advantage ReadyWide with mobility, floor lock & foot end controls
MC7WPML	Advantage ReadyWide Plus (650 lbs SWL) with mobility, floor lock & foot end controls
R9938	3 position molded pivot assist rail
R9939	Set of 3 position molded pivot assist rails
R9308H	Set of molded head end half side rails
R9308H-L	Left side molded head end half side rail
R9308H-R	Right side molded head end half side rail
R9518	Split side rail, 3/4 bed coverage
Q6749	Under bed lighting - factory installed
Q6701	Battery back up (Standard model) - factory installed
Q6770	Battery back up (Plus model) - factory installed
Q6740	Floor placement arm
Q6741	Retractable power cord cable reel & placement arm
Q6690	Bed extender to 84"
Q6730	Patient helper, trapeze bar, chain and hand grip
Q6731	Trapeze bed mounting bracket
Q6703	Bed and deck angle indicator
Q6699	Bed transportation system
Q6531	Honeycomb roller bumper kit

Contact your local Span Representative
or call Beds/Furniture Customer Service:
800-268-8671

ADVANTAGE
READYWIDE SERIES

Span's low-bed design with advanced clinical positioning and ReadyWide® integrated expansion

Also Available




Auto Contour

ReadyWide® Integrated Expansion



← Tool-less, "two-click" expansion. →

Advanced Clinical Positioning

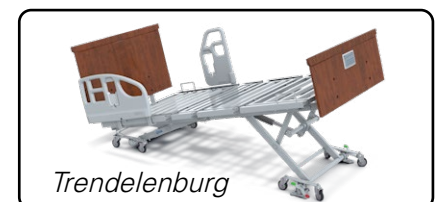

Trendelenburg




Reverse Trendelenburg

ADVANTAGE

READYWIDE SERIES

The Advantage ReadyWide Series with advanced features such as mobility at any height, easy optional sleep surface expansion to 42" (1067 mm), a multi-function carrier system and a safe working load of 500 lbs / 227 kg (Standard model) or 650 lbs / 295 kg (Plus model).

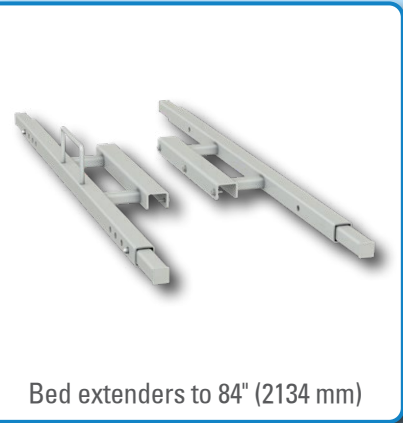

Bed extenders to 84" (2134 mm)



Trapezoidal slat profile provides unmatched strength, ventilation, and ease of cleaning.



Simple, push-button ReadyWide® operation: Click once to pull expanders out to 39" (991 mm) wide. Click again for 42" (1067 mm) width.



Optional, moulded quick release plastic half-head assist rails.



Optional 3-position side assist rail offers ingress-egress support.



In addition to the hand pendant controller, a foot-end control with staff-only lockout is standard.



Comfort Chair



Optional bed carrier facilitates easy under-bed service and bed passage on its side through narrow doorways.



Optional full 360° mobility and lockability at any height for safer, ergonomic bed relocation.



Optional dual-position, dual-intensity underbed lighting.

Smart Stop® automatically stops lowering when deck reaches 12.5" (318 mm) from floor.

Helps ensure safer, more ergonomically correct transfers. Positions typical user at ideal height* and leg angle for egress/regress.
(*12.5" / 318 mm plus height of mattress)

