
# Geo-Mattress Pro® RP

## Pressure Redistribution Raised Perimeter Mattress

### Affordability and added safety

Single layer, zoned construction delivers excellent value for budget-conscious facilities.

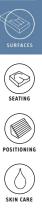

Supportive, raised perimeter helps prevent falls and entrapment by prompting user away from edge of bed.

ltem #	Description
PR8435RP-29	84"L x 35"W x 6"H (center), 8" (side)
PR8035RP-29	80"L x 35"W x 6"H (center), 8" (side)
PR7535RP-29	75"L x 35"W x 6"H (center), 8" (side)

attress 80"

- 5 years, non-prorated warranty
- Weight limit: 300 lbs.

Flammability Protection: Certified conformance to NFPA 101 (Life Safety Code); Cal. Tech Bulletin #117; 16 CFR 1632 and 1633; and BFD IX-11.

#### Primary Use:

• Prevention of Skin Breakdown

#### Patient Indications:


- General patient population
- Moderate High Risk for skin breakdown (e.g., Braden 13-26)
- Decreased mobility
- Repositions frequently by self or caregiver

#### Key Advantages:


- Geo-Matt® anti-shearing surface
- Heel Slope™
- Hidden zipper for infection control

#### Packed with value-added features:



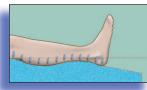



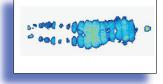

#### Segmented Geo-Matt® Design.

Clinically proven Geo-Matt<sup>®</sup> design creates a unique, anti-shearing top surface. 800 individually responsive cells in distinct head, torso, and foot zones enhance therapeutic support and comfort.

#### Exclusive Heel Slope<sup>®</sup>.

Subtle taper redistributes load to pressure-tolerant lower legs. Delivers heel pressures 27% lower than flat designs.

#### Fluid-proof, wipe-clean cover

with ultra-low moisture vapor transmission rate (MVTR). Bacteriostatically-treated top, durable tri-laminate bottom. Hidden zipper design for easier, more thorough cleaning and superior infection control.

#### Outstanding pressure redistribution

across the surface ensures low average of peak pressure at vulnerable sacrum and heels.

